

HPLC Application

ID No.: 10129

Aflatoxins

Column: Prodigy[™] 5 µm ODS-3 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4097-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Eluent C: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
	1	0	68	24	8

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: Fluorescence (FLUOR) @ 360/440 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy[™]

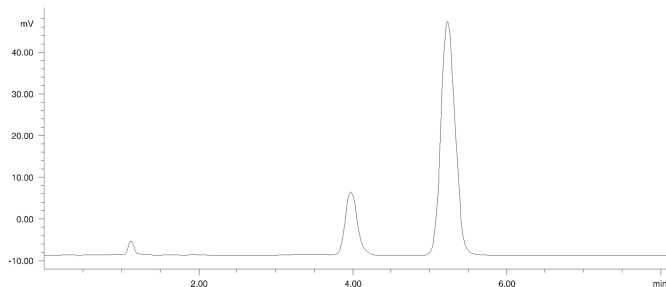

ANALYTES:

1 Aflatoxin M2

2 Aflatoxin M1

