

Fat-Soluble Vitamins

Column: Luna® 5 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	75	25

Flow Rate: 1.3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

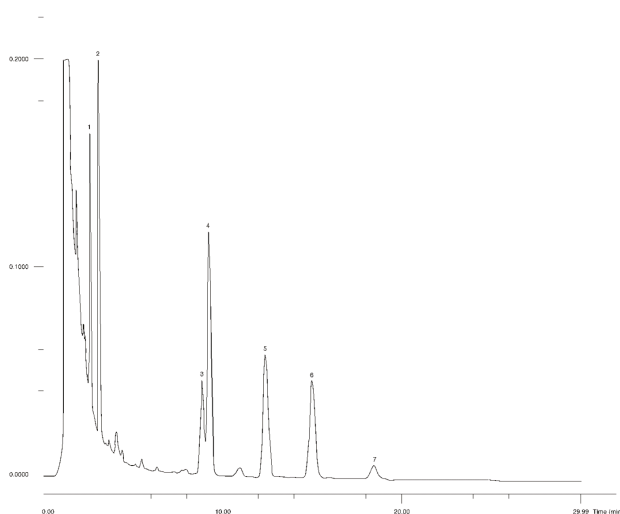

ANALYTES:

- 1 Vitamin A (trans-Retinol)
- 2 Vitamin A acetate
- 3 Vitamin D2 (Ergocalciferol)
- 4 Vitamin D3 (Cholecalciferol)
- 5 Vitamin E (alpha-Tocopherol)
- 6 Vitamin E acetate (alpha-Tocopherol acetate)
- 7 Vitamin K

