

Nucleotides

Column: Synergi™ 4 µm Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate pH 2.7

| Gradient Profile: | Step No. | Time (min) | Pct A |
|-------------------|----------|------------|-------|
| | 1 | 0 | 100 |

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282

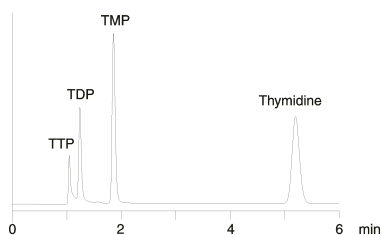

ANALYTES:

- 1 Thymidine triphosphate
- 2 Thymidine diphosphate
- 3 Thymidine monophosphate
- 4 Thymidine

