

HPLC Application

ID No.: 12952

Cyano Project - Competitor Comparisons - Corticosteroids

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4255-E0

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: 80:20 Methylene Chloride:Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305


- Holder Part No.: KJ0-4282



Products used in this application:



12952

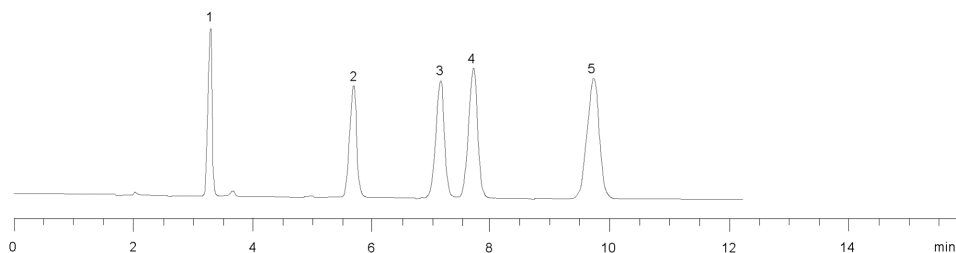

ANALYTES:

- 1 Desoxycorticosterone
- 2 Corticosterone
- 3 Cortisone
- 4 Prednisone
- 5 Prednisolone

