

HPLC Application

ID No.: 12953

Cyano Comparisons - Rhodamine Dyes

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4255-E0

Elution Type: Isocratic

Eluent A: 80:20 20 mM Potassium phosphate pH 3.0:Methanol

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305


- Holder Part No.: KJ0-4282



Products used in this application:



12953

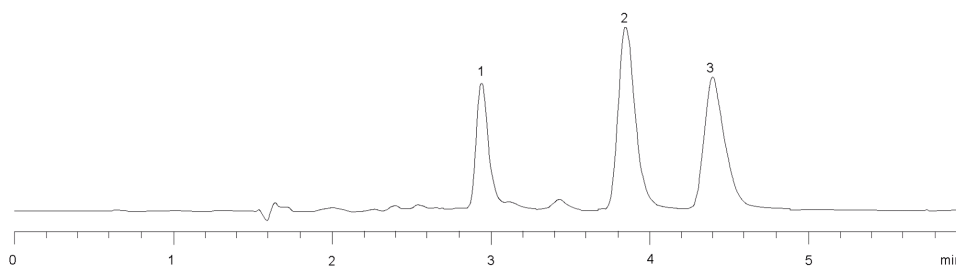

ANALYTES:

- 1 Rhodamine 123
- 2 Rhodamine B
- 3 Rhodamine 6G

