

HPLC Application

ID No.: 12958

Tricyclic Antidepressants on Luna CN

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4255-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 20 mM Potassium phosphate pH 4.0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282



Products used in this application:

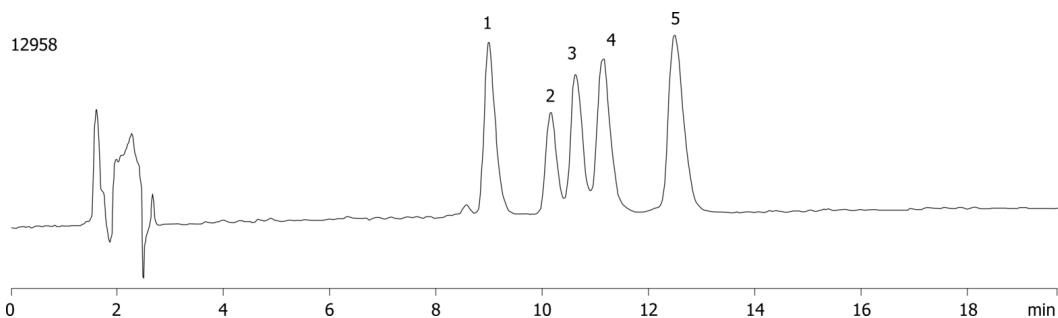

ANALYTES:

- 1 Doxepin
- 2 Desipramine
- 3 Nortriptyline
- 4 Amitriptyline
- 5 Protriptyline

