

GC Application

ID No.: 1328

Diesel Fuel Oil No. 2

Column: Zebron[™] ZB-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 30 meters x 0.25 mm x 0.25 µm

Order No: 7HG-G001-11

Oven Profile: 40°C (2 min) to 320°C @ 6°C/min (hold 5 min)


Carrier Gas: Constant Flow Hydrogen, 1.23 mL/min

Injection: Split 100:1 1 µL @ 250°C

Detection: Flame Ionization (FID) (325°C)



Products used in this application:

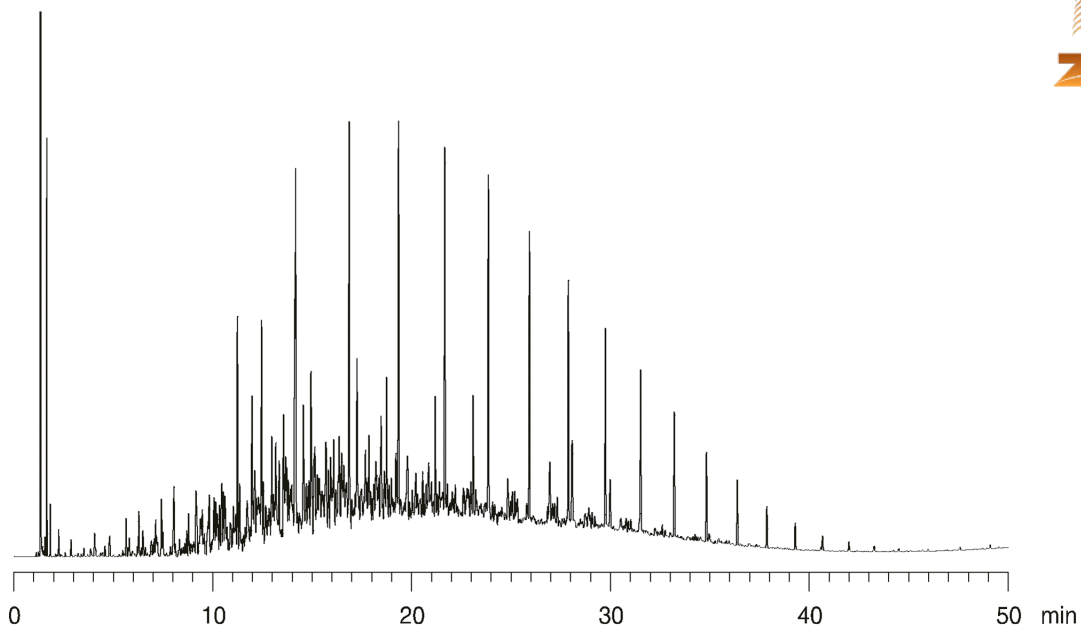

ANALYTES:

- 1 Diesel fuel (profile)

