

Phenolic compounds

Column: Prodigy™ 5 µm ODS-2 150 Å, LC Column 75 x 4.6 mm, Ea

Dimensions: 75 x 4.6 mm ID

Order No: 00C-3300-E0

Elution Type: Gradient

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	25	75
	2	11	30	70
	3	30	60	40
	4	30.1	100	0
	5	45	100	0

Flow Rate: 1 mL/min

Col. Temp.: 50 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 208 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

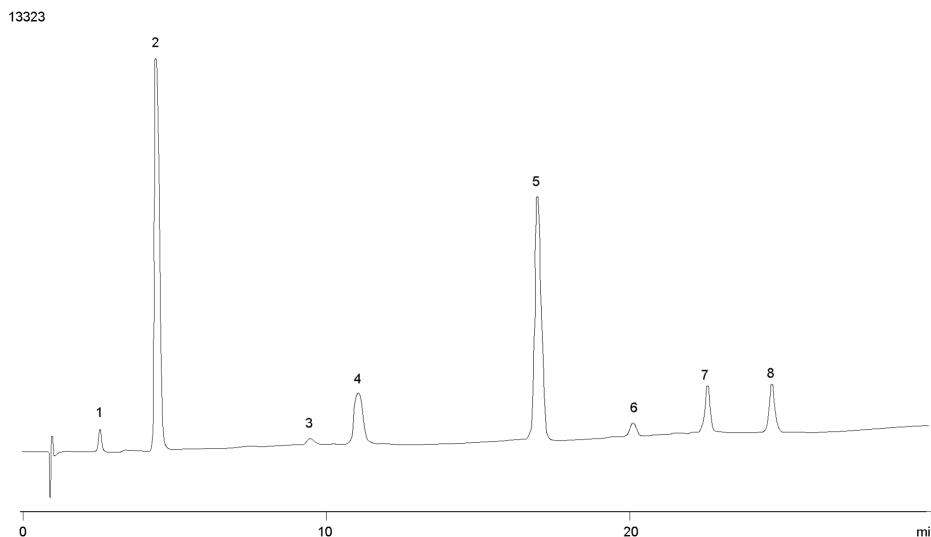

ANALYTES:

- 1 Phenol
- 2 o-Cresol
- 3 IPP
- 4 p,p-BPA
- 5 o,p-BPA
- 6 IPP2
- 7 BPX
- 8 COD

