

# HPLC Application

ID No.: 14049

## Asparagine

**Column:** Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-3126-E0

**Elution Type:** Isocratic



**Eluent A:** 1 mM Copper (II) sulfate in water

Gradient	Step No.	Time (min)	Pct A
<b>Profile:</b>	<b>1</b>	0	100

**Flow Rate:** 0.5 mL/min

**Col. Temp.:** 60 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

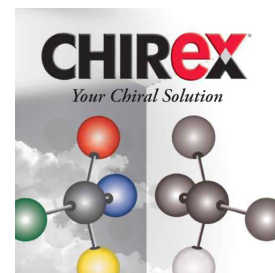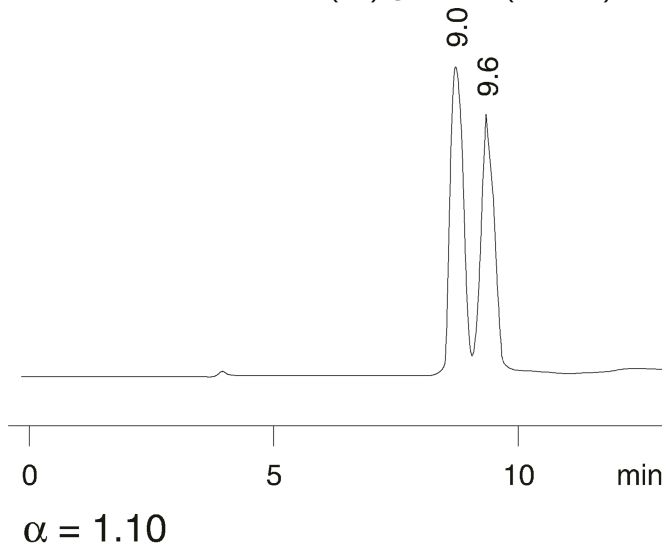

Products used in this application:

**CHIREX**®

### ANALYTES:

- 1** Asparagine  
Retention Time: 8.965 min
- 2** Asparagine  
Retention Time: 9.585 min