

HPLC Application

ID No.: 14279

Diquat by SCX

Column: Capcell Pak[®] 5 μ m SCX UG 80 \AA , LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: CH0-3882

Elution Type: Isocratic

Eluent A: 750mM KH₂PO₄ pH 2.7

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 308 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, SCX 4 x 3.0mm, 10/Pk Part No.: AJ0-4308

- Holder Part No.: KJ0-4282



Products used in this application:

Capcell Pak[®]

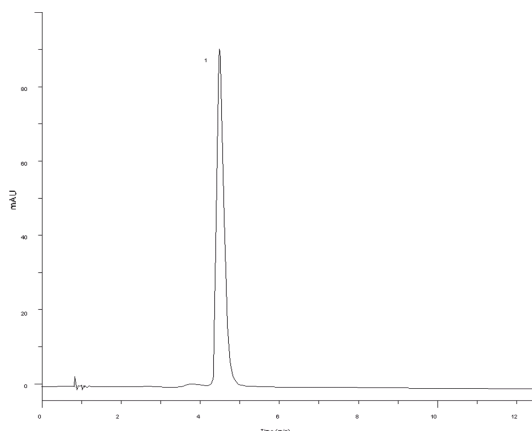

ANALYTES:

1 Diquat

