

GC Application

ID No.: 14577

Aroclors 1016 and 1260 on ZB5 15m x .32mm x .25um

Column: Zebron[™] ZB-5, GC Cap. Column 15 m x 0.32 mm x 0.25 μ m, Ea

Phase: 5% Phenyl 95% Dimethylpolysiloxane

Dimensions: 15 meters x 0.32 mm x 0.25 μ m

Order No: 7EM-G002-11


Oven Profile: 150°C to 300°C at 12°C/min

Carrier Gas: Constant Flow Helium, 1.7 mL/min

Injection: split 15:1:1 1 μ L @ 250°C

Detection: Electron Capture (ECD) (375°C)

14577

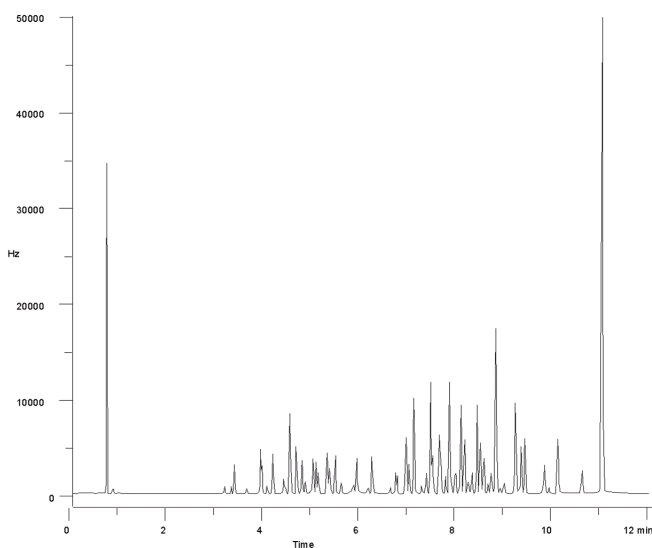

Products used in this application:



ANALYTES:

- 1 Aroclor 1016
- 2 Aroclor 1260