

HPLC Application

ID No.: 14766

siRNA on Jupiter Proteo

Column: Jupiter[®] 4 μ m Proteo 90 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4396-E0

Elution Type: Gradient

Eluent A: 100mM TEAA pH 7

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.5	96	4
	2	25	88	12

Flow Rate: 1 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282



Products used in this application:

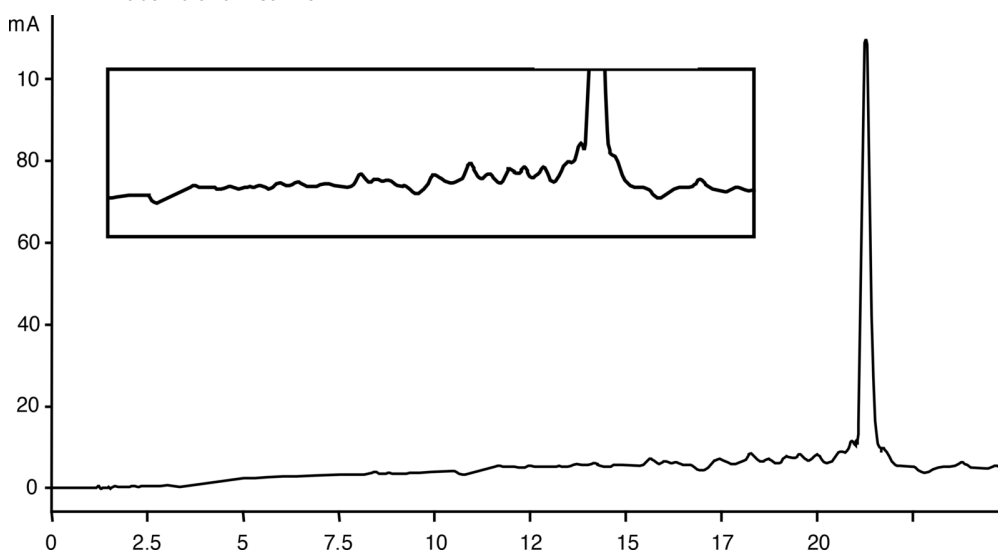

ANALYTES:

1 siRNA 19mer

