

## Aqueous Stability using Synergi Fusion-RP

**Column:** Synergi™ 4 µm Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4424-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM Potassium Phosphate buffer, pH 2.5

Gradient	Step No.	Time (min)	Pct A
<b>Profile:</b>	<b>1</b>	0	100

**Flow Rate:** 1 mL/min

**Col. Temp.:** 30 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



14891

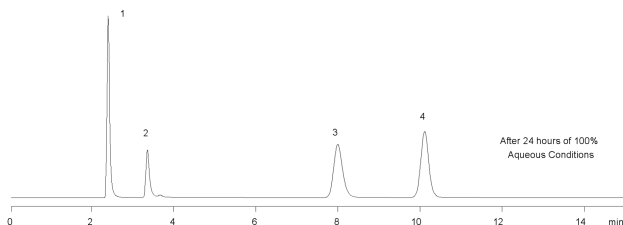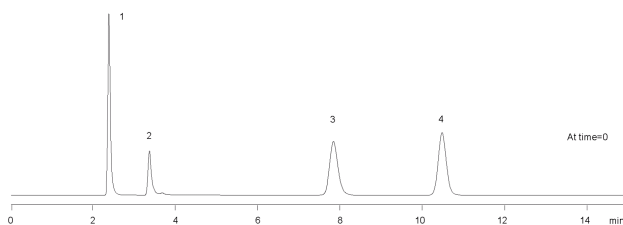

### ANALYTES:

- 1 Thiourea
- 2 Adenine
- 3 Guanosine-5-monophosphate
- 4 Thymine

