

## Caffeine & Metabolites, pH 4.5

**Column:** Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4435-E0

**Elution Type:** Gradient

**Eluent A:** 20mM Ammonium acetate, pH 4.5

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2
	2	15	80	20

**Flow Rate:** 1 mL/min

**Col. Temp.:** 25 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 273 nm (25 °C)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597

- Holder Part No.: KJ0-4282

15013



Products used in this application:

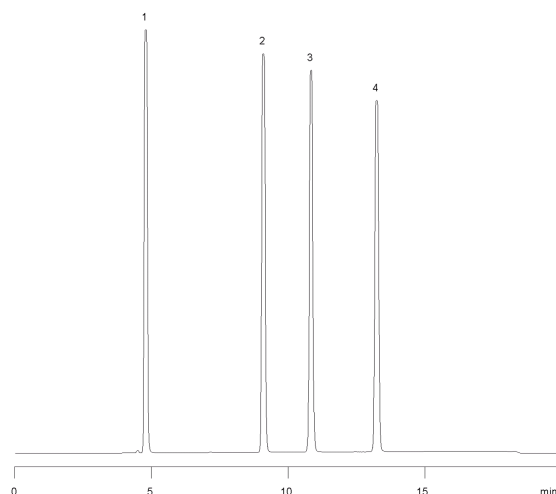

### ANALYTES:

- 1** Xanthine  
Retention Time: 4.791 min
- 2** Theobromine  
Retention Time: 9.092 min
- 3** Theophylline  
Retention Time: 10.833 min
- 4** Caffeine  
Retention Time: 13.235 min

