

HPLC Application

ID No.: 15021

Non-Polar Acids -- Analgesics

Column: Gemini[®] 5 μ m C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 20mM Phosphate Buffer, pH2.5

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597

- Holder Part No.: KJ0-4282

15021

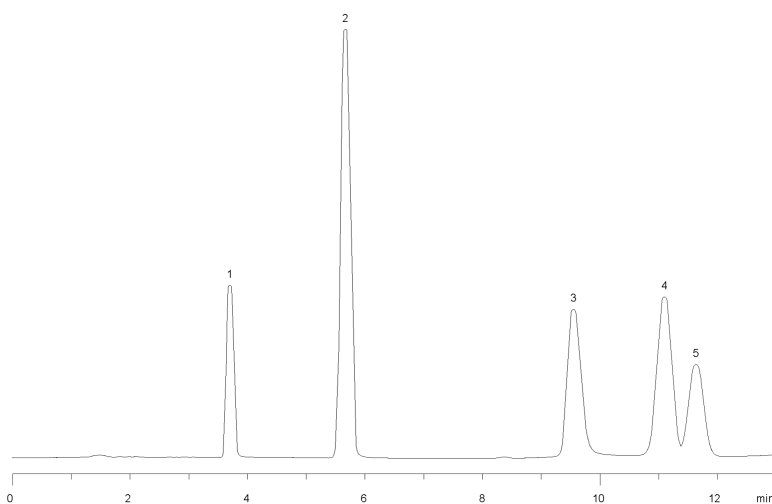

ANALYTES:

- 1 Ethyl paraben
- 2 Naproxen
- 3 Diflunisal
- 4 Indomethacin
- 5 Ibuprofen



Products used in this application:

