

GC Application

ID No.: 15426

Butane Purity ZB1 60x0.32x0.25

Column: Zebtron™ ZB-1, GC Cap. Column 60 m x 0.32 mm x 0.25 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 60 meters x 0.32 mm x 0.25 µm

Order No: 7KM-G001-11

Oven Profile: 40°C Isothermal

Carrier Gas: Constant Flow Hydrogen, 6.5 mL/min


Injection: split 10:1 5 µL @ 200°C

Detection: Flame Ionization (FID) (200°C)

15426



Products used in this application:

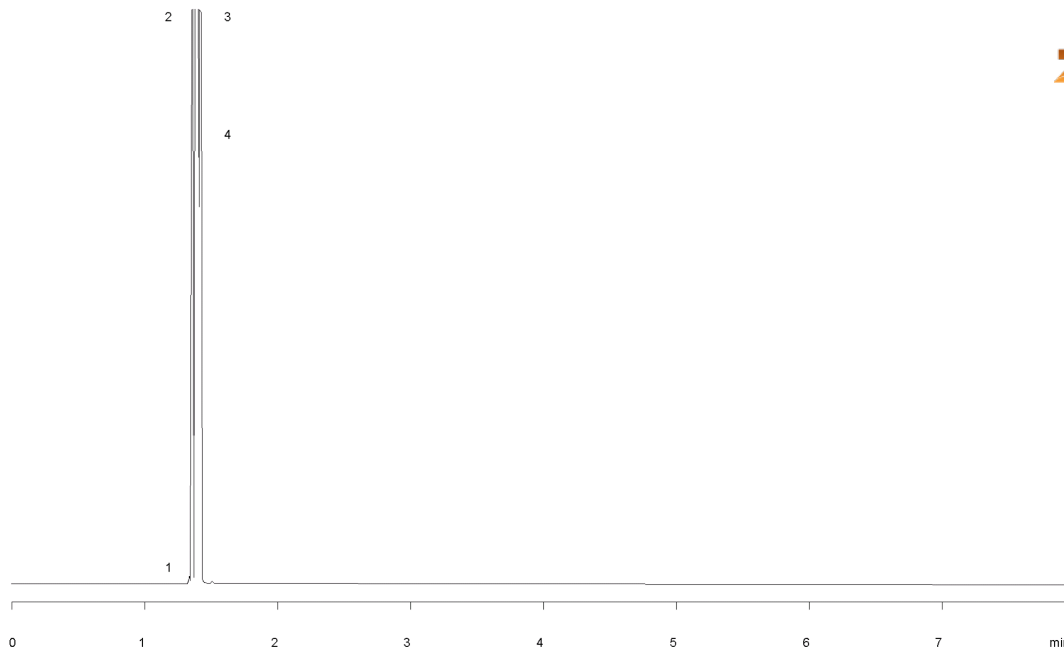

ANALYTES:

- 1 Ethane
- 2 Propane
- 3 Isobutane
- 4 Butane

