


Organo Phosphorous Pesticides Mix 1 on ZB-5MS

Column: Zebron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 100°C to 140°C at 20°C/min to 200°C at 5°C/min to 300°C at 25°C/min for 5 min
Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: split 20:1 1 µL @ 250°C
Detection: Mass Selective (MSD) (240°C)

15435



Products used in this application:

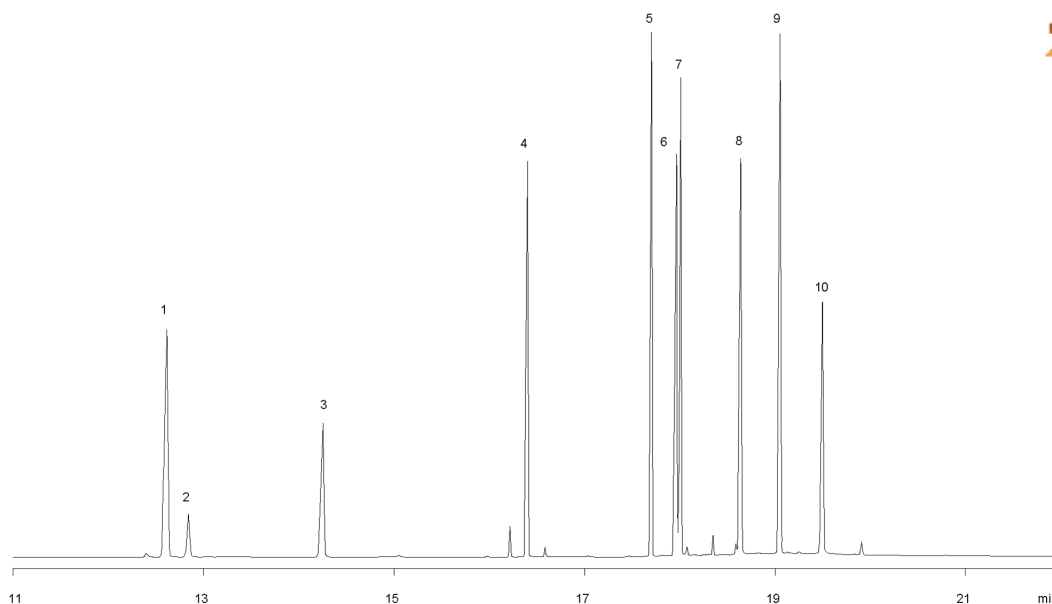

Organo Phosphorous Pesticides Mix 1 on ZB-5MS

ANALYTES:

- 1 Terbufos (Terbuphos)
- 2 Phosphamidon (isomer)
- 3 Phosphamidon
- 4 Chlorfenvinphos
- 5 Ethion
- 6 Famphur (Famophos)
- 7 Carbophenothion
- 8 Phosmet
- 9 Leptophos
- 10 Azinphos-ethyl

