

# HPLC Application

ID No.: 15741

## Acyclovir and Ganciclovir analysis using Synergi Fusion-RP

**Column:** Synergi<sup>™</sup> 4  $\mu$ m Fusion-RP 80  $\text{\AA}$ , LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4424-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM Ammonium Formate, pH 3.5

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



15741

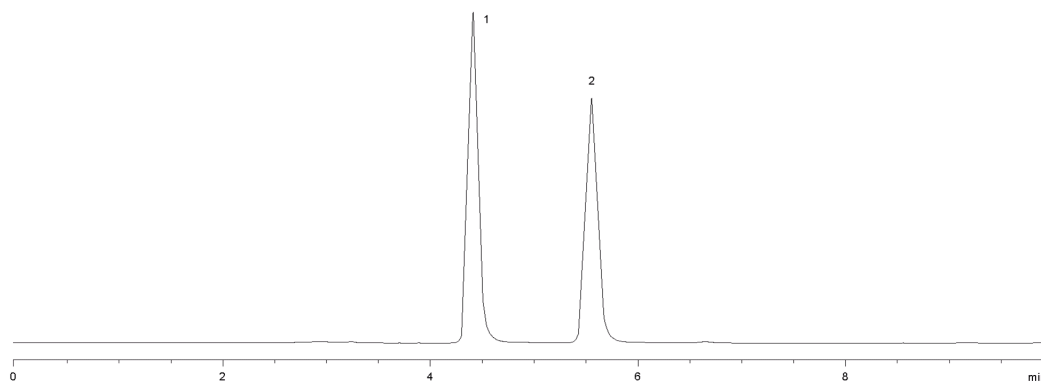

### ANALYTES:

- 1 Acyclovir
- 2 Ganciclovir

