

Massachusetts EPH - Aliphatic Fraction on ZB-5MS

Column: Zebtron™ ZB-5MS, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea

Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane

Dimensions: 30 meters x 0.32 mm x 0.25 µm

Order No: 7HM-G010-11

Oven Profile: 60°C for 1 min to 290°C at 8°C/min for 6.75 min

Carrier Gas: Constant Flow Helium, 3 mL/min


Injection: Splitless :1 2 µL @ 285°C

Detection: Flame Ionization (FID) (315°C)

Analyst Note: Analytes are 50ppm and IS and Surrogates are 40ppm for fractionation.



Products used in this application:

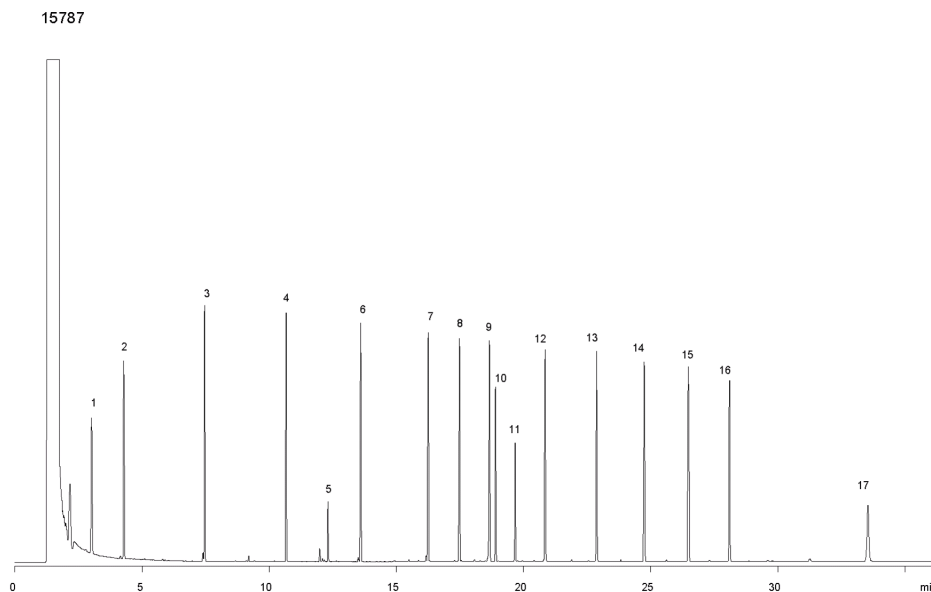

Massachusetts EPH - Aliphatic Fraction on ZB-5MS

ANALYTES:

- 1 C9
- 2 C10
- 3 C12
- 4 C14
- 5 Butylhydroxytoluene
- 6 C16
- 7 C18
- 8 C19
- 9 C20
- 10 5-alpha-Androstane (internal standard)
- 11 1-Chloro-octadecane (surrogate)
- 12 C22
- 13 C24
- 14 C26
- 15 C28
- 16 C30
- 17 C36

