

Water Reproducibility of ZB-WAX Plus 30x0.25x0.25

Column: Zebtron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea

Phase: Polyethylene Glycol (PEG)

Dimensions: 30 meters x 0.25 mm x 0.25 μ m

Order No: 7HG-G013-11

Oven Profile: 35°C for 5 min to 85°C at 10°C/min to 200°C at 25°C/min for 1 min

Carrier Gas: Constant Flow Helium, 1.4 mL/min


Injection: Split 30:1 0.2 μ L @ 140°C

Detection: Flame Ionization (FID) (200°C)

15817



Products used in this application:

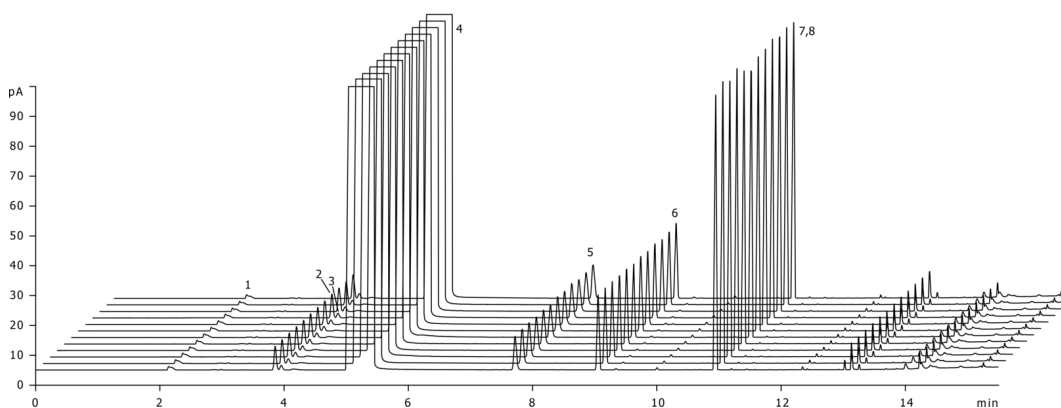

ANALYTES:

- 1 Acetaldehyde
- 2 Ethyl acetate
- 3 Methanol
- 4 Ethanol
- 5 Propanol
- 6 Isobutanol
- 7 2-Methylbutanol
- 8 3-Methylbutanol

