

Jupiter® 300 C4 - Comparison of Mobile Phase Conditions

Column: Jupiter® 5 µm C4 300 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4167-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 90% ACN 0.09% TFA in Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	25	45	55


Flow Rate: 1 mL/min

Col. Temp.: 45 °C


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

Analyst Note: Reversed phase runs of PEGylated carbonic anhydrase and PEGylated ribonuclease A using different organic mobile phase compositions. In other applications, the addition of isopropanol does not improve resolution of different PEGylated forms.

16194

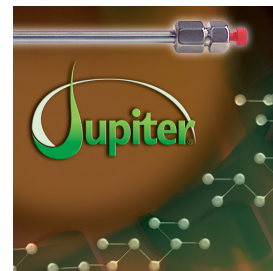

Products used in this application:

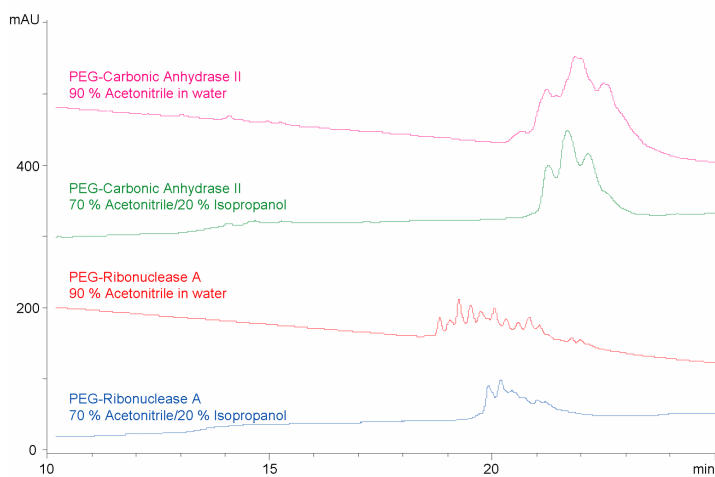

ANALYTES:

- 1 PEGylated Carbonic anhydrase or PEGylated Ribonuclease A

