


ASTM Method D2887-97: Boiling Range Distribution of Petroleum Fractions on ZB-MR2

Column: Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μm, Ea
Phase: Unknown or Proprietary
Dimensions: 30 meters x 0.25 mm x 0.2 μm
Order No: 7HG-G017-10
Oven Profile: 70 °C to 320 °C at 8 °C/min.

Carrier Gas: Constant Flow Helium, 4 mL/min
Injection: Split 10:1 0.5 μL @ 310°C
Detection: Flame Ionization (FID) (320°C)
Analyst Note: Analytes 100ppm in Hexane



Products used in this application:

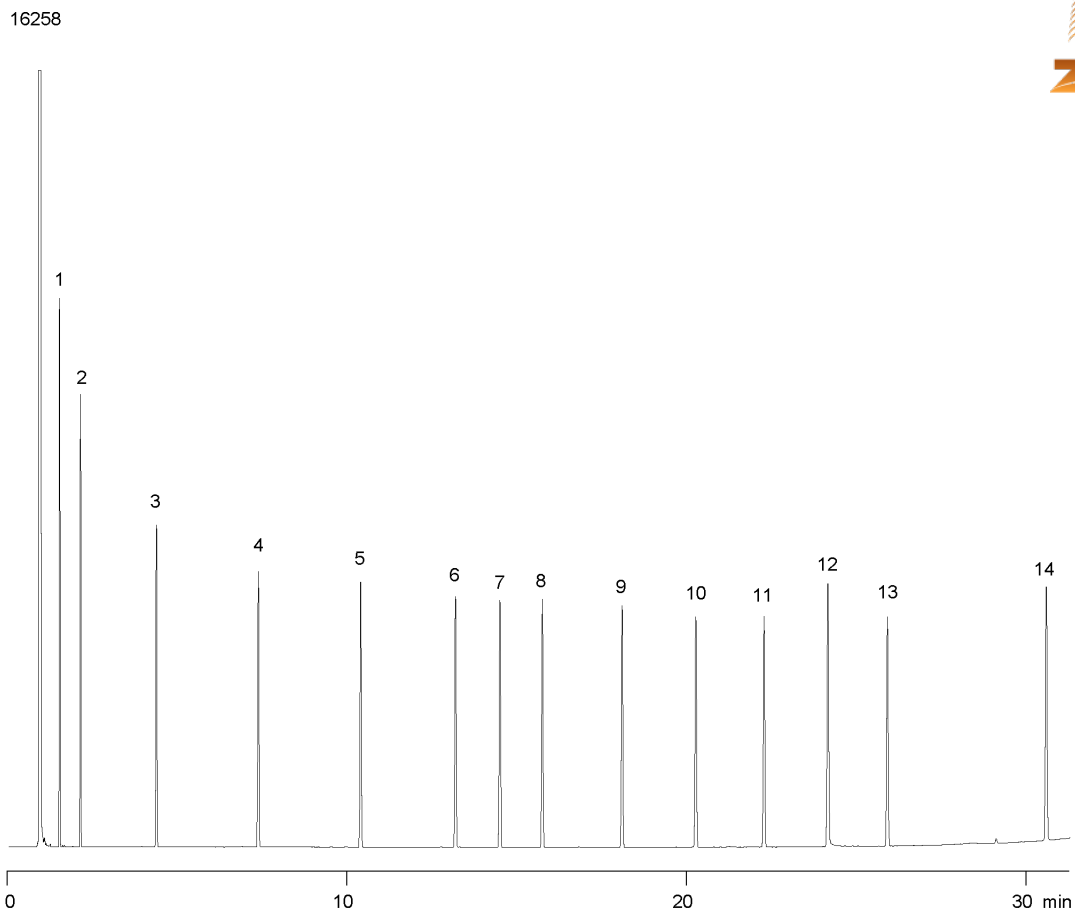

ASTM Method D2887-97: Boiling Range Distribution of Petroleum Fractions on ZB-MR2

ANALYTES:

- 1 n-Nonane (C9)
- 2 n-Decane (C10)
- 3 n-Dodecane (C12)
- 4 n-Tetradecane (C14)
- 5 n-Hexadecane (C16)
- 6 n-Octadecane (C18)
- 7 n-Nonadecane (C19)
- 8 n-Eicosane (C20)
- 9 n-Docosane (C22)
- 10 n-Tetracosane (C24)
- 11 n-Hexacosane (C26)
- 12 n-Octacosane (C28)
- 13 n-Tricosane (C30)
- 14 n-Hexatricosane (C36)

