

HPLC Application

ID No.: 16399

Luna HILIC Test Mix (AL0-8317)

Column: Luna® 5 µm HILIC 200 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4450-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 100mM Ammonium formate, pH 3.2

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, HILIC 4 x 3.0 mm, 10/Pk Part No.: AJ0-8329

- Holder Part No.: KJ0-4282



Products used in this application:

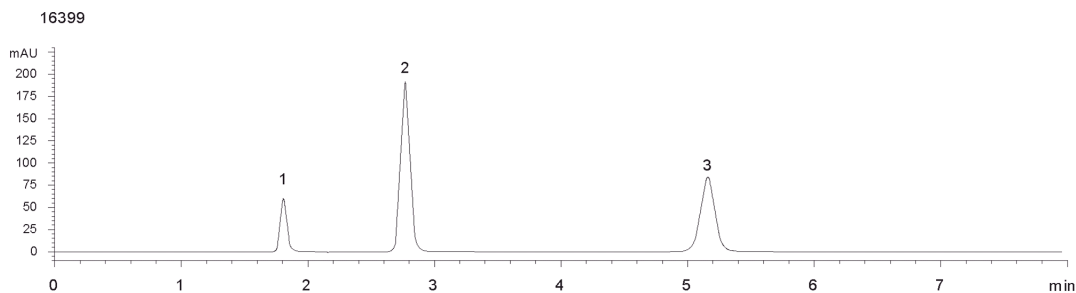

ANALYTES:

- 1 Toluene (void volume t₀)
- 2 Uracil
- 3 Cytosine

