

GC Application

ID No.: 16472

Aroclor 1246 on ZB-XLB

Column: Zebtron[™] ZB-XLB, GC Cap. Column 15 m x 0.25 mm x 0.25 µm, Ea

Phase: Proprietary XLB Phase

Dimensions: 15 meters x 0.25 mm x 0.25 µm

Order No: 7EG-G019-11

Oven Profile: 150 °C to 300 °C @ 15 °C/min.

Carrier Gas: Constant Flow Hydrogen, 2 mL/min

Injection: Split 10:1 1 µL @ 240°C


Detection: Electron Capture (ECD) (320°C)

Analyst Note: Aroclor 1246 concentration 1 ug/ml

App ID 16472



Products used in this application:

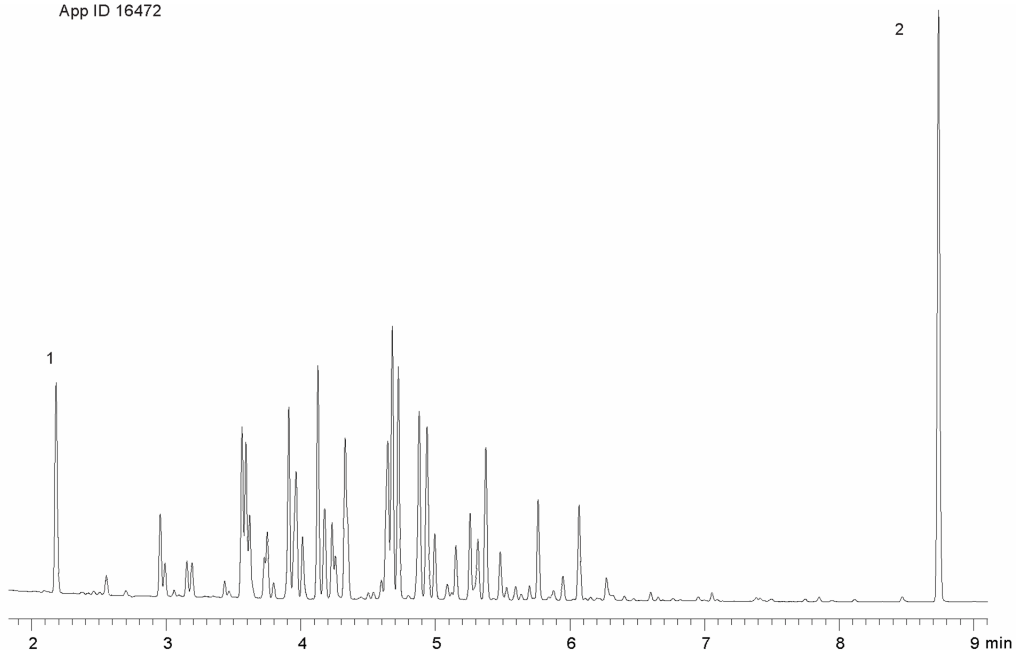

ANALYTES:

- 1 Tetrachloro-m-xylene (internal standard)
- 2 Decachlorobiphenyl (internal standard)

