

GC Application

ID No.: 17569

NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix on ZB-MR2

Column: Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μm, Ea

Phase: Proprietary Pesticides Phase

Dimensions: 30 meters x 0.25 mm x 0.2 μm

Order No: 7HG-G017-10

Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min

Carrier Gas: Constant Flow Helium, 0.8 mL/min

Injection: On-Column :1 0.5 μL @ 43°C

Detection: Mass Selective (MSD) (335°C)


Analyst Note: Analytes 10 ppm in Acetone



Products used in this application:



17569

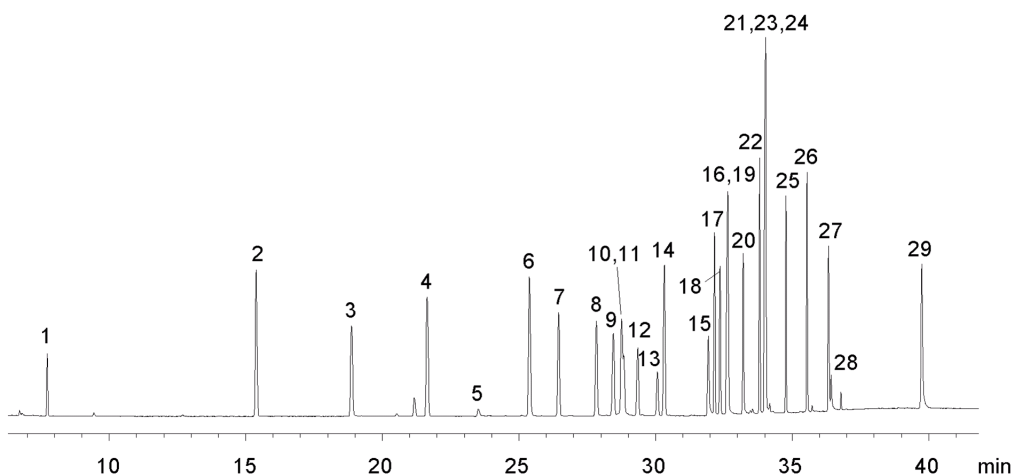

NSO (Nitrogen/Sulfur/Oxygen) Pesticides Mix on ZB-MR2

ANALYTES:

- 1 Metaldehyde
- 2 Biphenyl
- 3 o-Phenylphenol
- 4 Diphenylamine
- 5 Ethoxyquin
- 6 Pyrimethanil
- 7 Pirimicarb
- 8 Metalaxyl
- 9 Carbaryl
- 10 Methiocarb
- 11 Diethofencarb
- 12 MGK-264 (isomer 1)
- 13 MGK-264 (isomer 2)
- 14 Cyprodinil
- 15 Thiabendazole
- 16 Isoprothiolane
- 17 Buprofezin
- 18 Bupirimate
- 19 Flusilazole
- 20 Diniconazole
- 21 Oxadixyl
- 22 Benalaxyl
- 23 Propargite
- 24 Piperonyl butoxide
- 25 Fenpropathrin (Danitol)
- 26 Pyriproxifen
- 27 Bitertanol (isomer 1)
- 28 Bitertanol (isomer 2)
- 29 Azoxystrobin

