

GC Application

ID No.: 17611

Propoxyphene by MSD on a ZB-XLB

Column: Zebron[™] ZB-XLB, GC Cap. Column 10 m x 0.18 mm x 0.18 μm, Ea

Phase: Proprietary XLB Phase

Dimensions: 10 meters x 0.18 mm x 0.18 μm

Order No: 7CD-G019-08

Oven Profile: 170 °C to 320 °C @ 20 °C/min

Carrier Gas: Constant Flow Helium, 55 cm/sec

Injection: Split 10:1 0.1 μL @ 280°C

Detection: Mass Selective (MSD) (320°C)


Analyst Note: Analytes are 24ppm
Sim ions:
(Quant ion*) D11-Propoxyphene (*64,198,213): Propoxyphene (*58,193,208)



Products used in this application:



17611

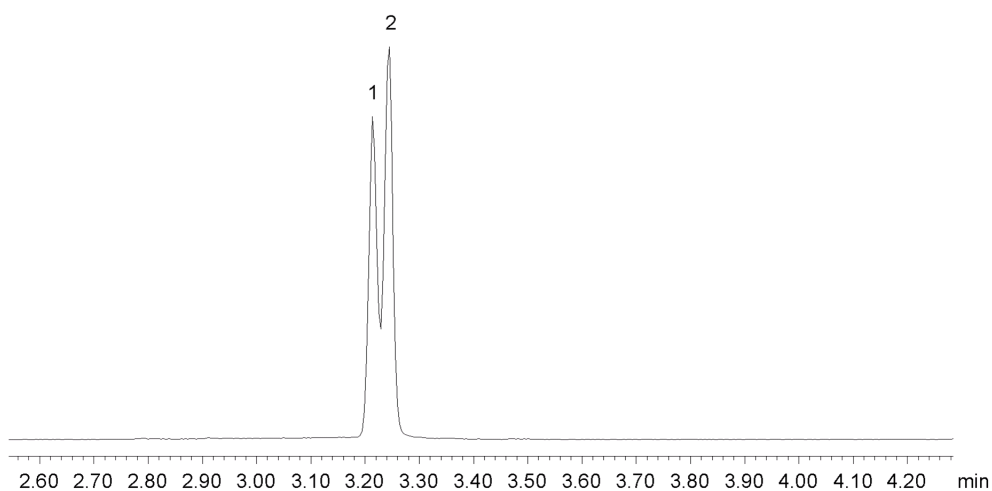

ANALYTES:

- 1 Propoxyphene-d11
- 2 Propoxyphene

