

## Expanded amphetamines by MSD on a ZB-1MS

**Column:** Zebron<sup>™</sup> ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea

**Phase:** 100% Dimethylpolysiloxane

**Dimensions:** 10 meters x 0.18 mm x 0.18 µm

**Order No:** 7CD-G011-08

**Oven Profile:** 100 °C to 320 °C @ 25 °C/min

**Carrier Gas:** Constant Flow Helium, 1 mL/min

**Injection:** Split 50:1 1 µL @ 280°C


**Detection:** Mass Selective (MSD) (230°C)

**Analyst Note:** Analytes are 40ppm derivatized with HFBA

Sim ions:

(Quant ion\*) D11-Amphetamine (\*244,128): Amphetamine (\*240,118,117): D14-Methamphetamine (\*261,213):  
Methamphetamine (\*254,210,118): D5-MDA (\*167,380,244): MDA (\*162,375,240): D5-MDMA (\*258,213,164):

17618

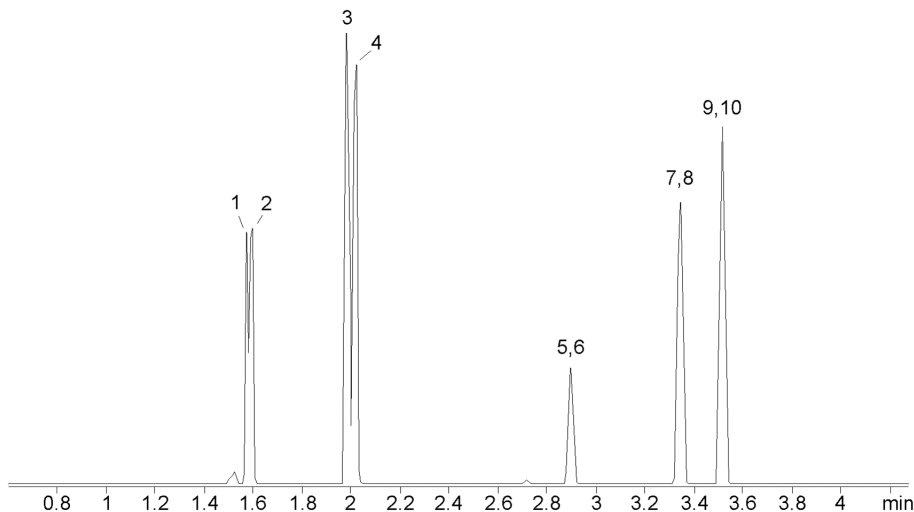

Products used in this application:



## Expanded amphetamines by MSD on a ZB-1MS

### ANALYTES:

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine
- 5 MDA-d5
- 6 MDA
- 7 MDMA-d5
- 8 MDMA
- 9 MDEA-d5
- 10 MDEA

