

GC Application

ID No.: 17805

Benzoyllecgonine by MSD on a ZB-Drug1

Column: Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea

Phase: Proprietary Drug Phase

Dimensions: 10 meters x 0.18 mm x 0.18 µm

Order No: 7CD-G023-08

Oven Profile: 235°C to 320°C @ 25 °C/min

Carrier Gas: Constant Flow Helium, 55 cm/sec

Injection: Split 10:1 1 µL @ 260°C

Detection: Mass Selective (MSD) (230°C)


Analyst Note: Analytes are 24ppm derivatized with BSTFA
Sim ions:
(Quant ion*) D3-Benzoyllecgonine (*243,259,364): Benzoyllecgonine (*240,256,361)



Products used in this application:



17805

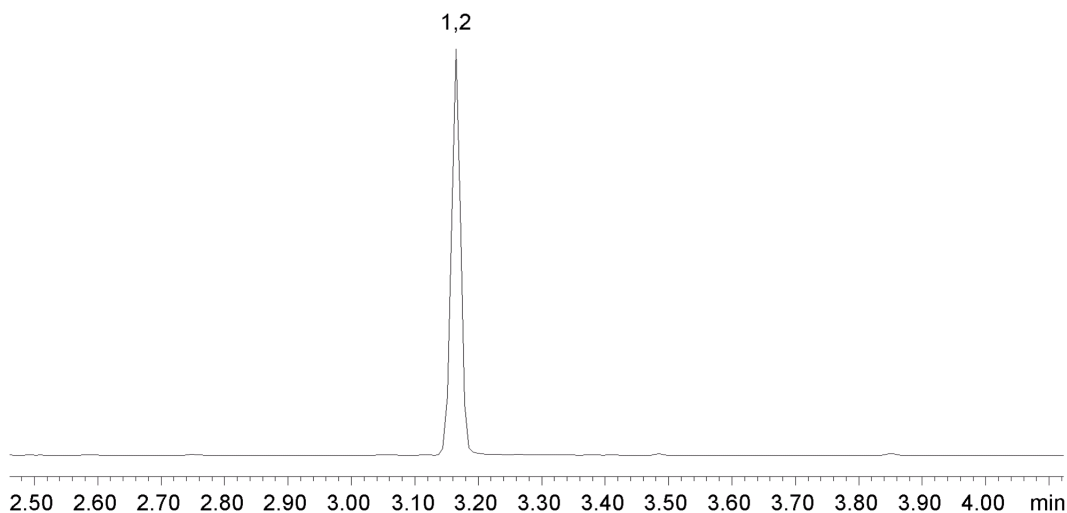

ANALYTES:

- 1 Benzoyllecgonine-d3
- 2 Benzoyllecgonine