

HPLC Application

ID No.: 18102

Brompheniramine on Lux 5 µm Amylose-2 250 x 4.6mm in NP

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: 0.1 % DEA in Hexane

Eluent B: 0.1% DEA in IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	10	90	10

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

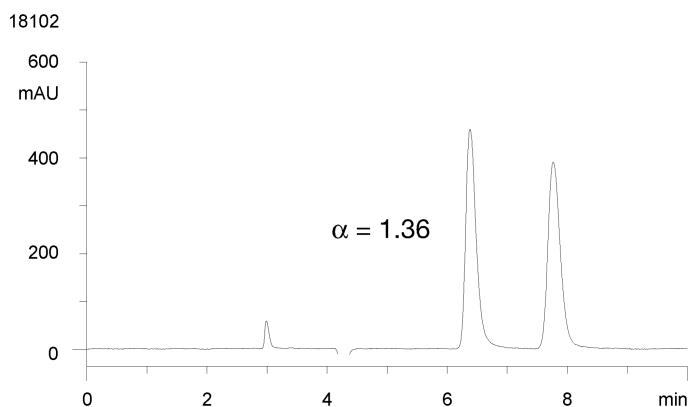

ANALYTES:

- 1** Brompheniramine
Retention Time: 6.38 min
- 2** Brompheniramine
Retention Time: 7.77 min

