

Acebutolol on Lux 5 µm Amylose-2 250 x 4.6mm in PO-95%CH3CN


Column: N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1 % DEA in CH3CN**Eluent B:** 0.1% DEA in IPA**Gradient**

Profile:	Step No.	Time (min)	Pct A	Pct B
	1	10	95	5

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

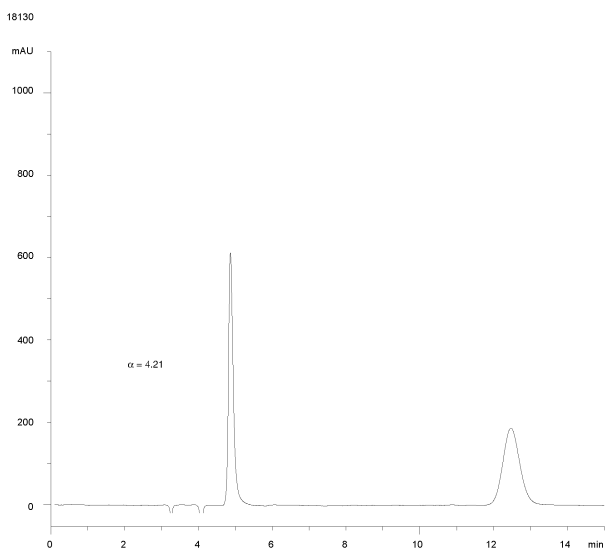

ANALYTES:

1 Acebutolol

Retention Time: 4.87 min

2 Acebutolol

Retention Time: 12.47 min

