

Column Test Mix on Kinetex 2.6µm C18 150 x 4.6mm ID

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Isocratic

Eluent A: water

Eluent B: acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	30	70

Flow Rate: 1.8 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.


- SecurityGuard ULTRA Cartridges, UHPLC C18 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8768

- Holder Part No.: AJ0-9000

18185



Products used in this application:

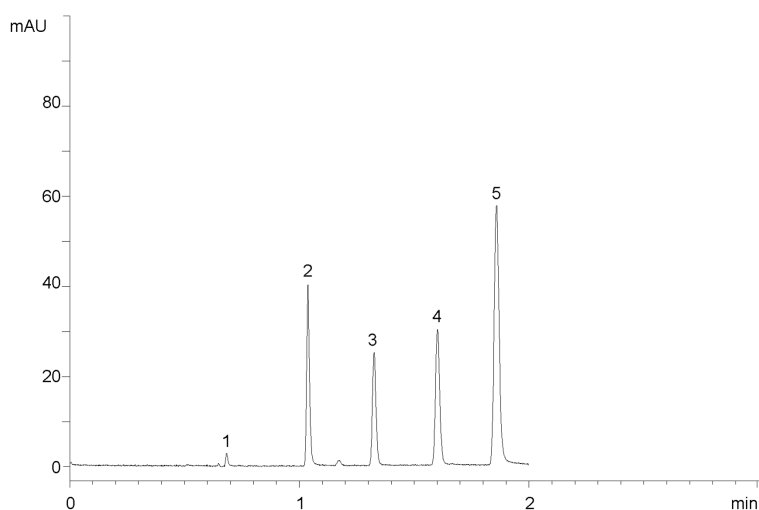

ANALYTES:

- 1** Uracil
Retention Time: 0.684 min
- 2** Acetophenone
Retention Time: 1.038 min
- 3** Benzene
Retention Time: 1.325 min
- 4** Toluene
Retention Time: 1.601 min
- 5** Naphthalene
Retention Time: 1.858 min

