

Ketones on Kinetex 2.6µm C18 150 x 4.6mm ID (1.7 mL/min)

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Gradient

Eluent A: water

Eluent B: acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1.21	95	5
	3	12.99	5	95
	4	15	5	95
	5	15.01	95	5
	6	20.13	95	5

Flow Rate: 1.7 mL/min

Col. Temp.: 45 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 258 nm (25 °C)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C18 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8768


- Holder Part No.: AJ0-9000



Products used in this application:



18278

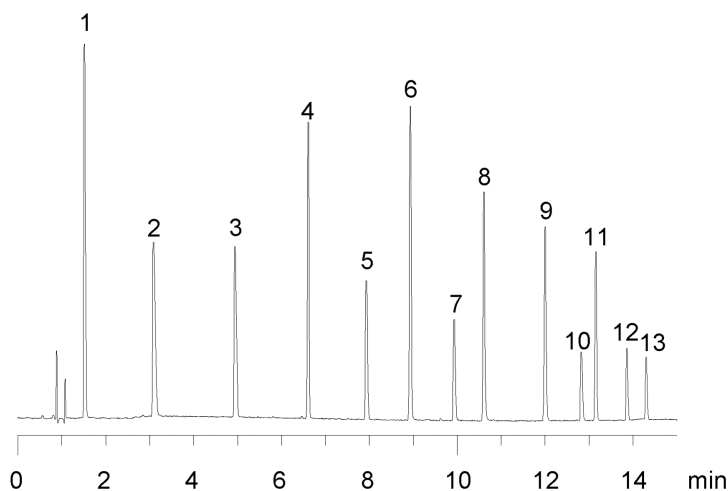

Ketones on Kinetex 2.6 μ m C18 150 x 4.6mm ID (1.7 mL/min)

ANALYTES:

- 1 Acetone
- 2 2-Butanone (MEK)
- 3 2-Pentanone
- 4 Acetophenone
- 5 2-Heptanone
- 6 Butyrophenone
- 7 2-Nonanone
- 8 Hexanophenone
- 9 Octanophenone
- 10 2-Tridecanone
- 11 Decanophenone
- 12 2-Pentadecanone
- 13 2-Hexadecanone

