

Triazine Pesticides (EPA 619) Mix # 1 on ZB-5MSi


Column: Zebtron™ ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μm
Order No: 7HG-G018-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 μL @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone



Products used in this application:



18328

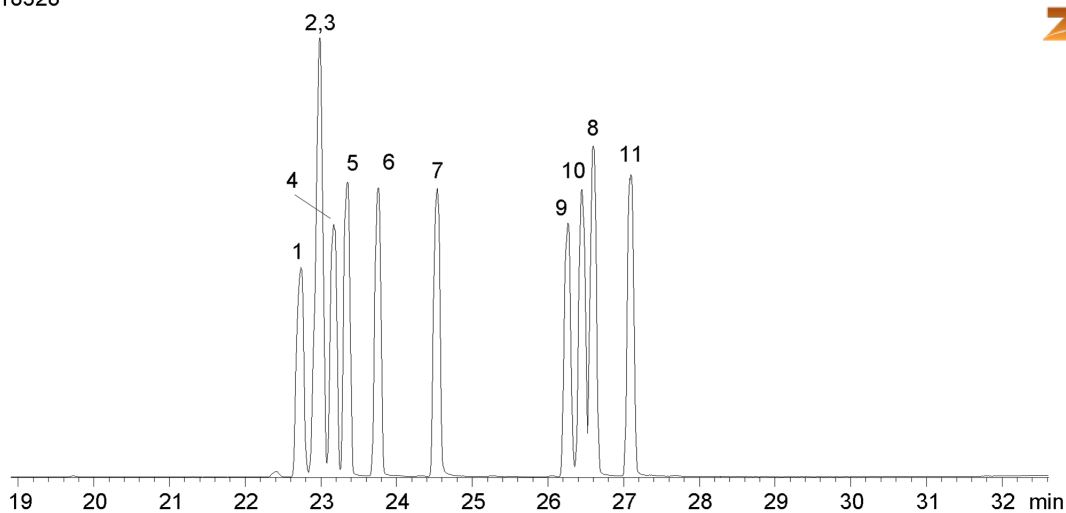

Triazine Pesticides (EPA 619) Mix # 1 on ZB-5MSi

ANALYTES:

- 1 Atraton
- 2 Prometon
- 3 Simazine
- 4 Atrazine
- 5 Propazine
- 6 Terbutylazine (Terbutylazine)
- 7 Secbumetone
- 8 Prometryn
- 9 Simetryn
- 10 Ametryn
- 11 Terbutryn

