

## Organophosphorous Pesticides (EPA 8141) Mix # 4 on ZB-35


**Column:** Zebtron™ ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** 35% Phenyl 65% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G003-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min  
**Carrier Gas:** Constant Flow Helium, 0.8 mL/min  
**Injection:** On-Column :1 0.5 µL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes 10 ppm in Acetone



Products used in this application:



18395

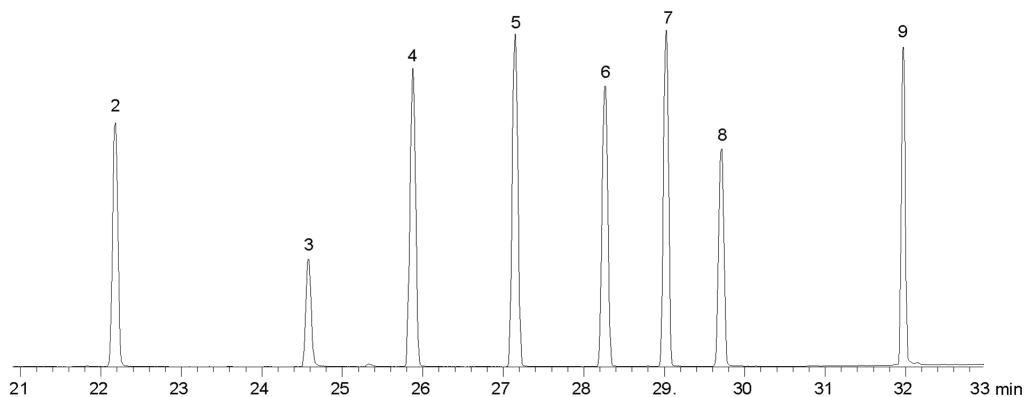

### ANALYTES:

- 1 Trichlorfon
- 2 Thionazin (Zinophos)
- 3 Dicrotophos
- 4 Fonofos
- 5 Dichlofenthion
- 6 Chlorpyrifos-methyl
- 7 Aspon
- 8 Fenitrothion
- 9 Crotoxyphos

