


GC Application

ID No.: 18407

Nitrogen and Phosphorous Pesticides (EPA 507) Mix # 3 on ZB-35

Column: Zebron[™] ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 35% Phenyl 65% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G003-11
Oven Profile: 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: On-Column :1 0.5 μ L @ 43°C
Detection: Mass Selective (MSD) (335°C)
Analyst Note: Analytes 10 ppm in Acetone

18407

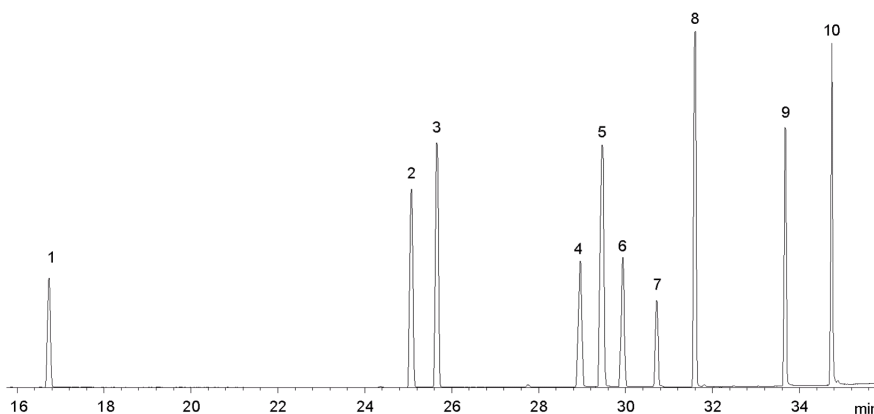

Products used in this application:



Nitrogen and Phosphorous Pesticides (EPA 507) Mix # 3 on ZB-35

ANALYTES:

- 1 Vernolate
- 2 Terbufos (Terbuphos)
- 3 Diazinon
- 4 Metribuzin
- 5 Metolachlor
- 6 MGK-264 (isomer 1)
- 7 MGK-264 (isomer 2)
- 8 Butachlor
- 9 Carboxin
- 10 Norflurazon

