

6-MAM, PCP and Methadone on Kinetex 2.6µm C18 50 x 2.1mm ID with Acquity UPLC MS/MS (1 mL/min)

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% Formic acid in H2O

Eluent B: 0.1% Formic acid in ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	1	5	95
	3	1.4	5	95
	4	1.41	5	95
	5	2	90	10

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: Mass Spectrometer (MS) @ amu (110 °C)


Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C18 for 2.1mm ID Columns, 3/Pk Part No.: AJ0-8782

- Holder Part No.: AJ0-9000



Products used in this application:

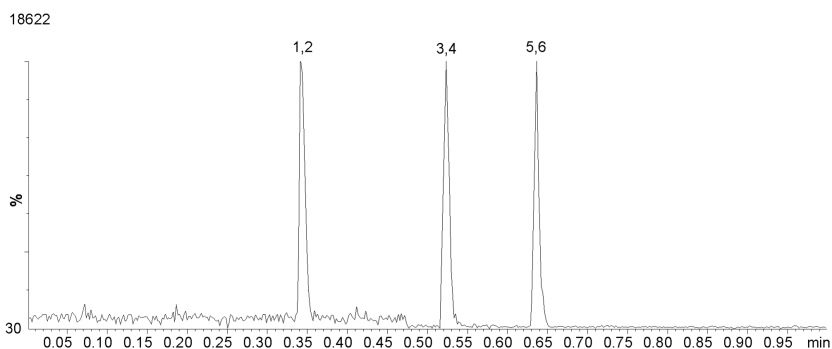

ANALYTES:

- 1 6-MAM
Retention Time: 0.35 min
- 2 6-MAM-d3
Retention Time: 0.35 min
- 3 PCP
Retention Time: 0.52 min
- 4 PCP-d5
Retention Time: 0.52 min
- 5 Methadone
Retention Time: 0.64 min
- 6 Methadone-d9
Retention Time: 0.64 min

