


FMOC-DL-Proline-(boc)-OH on Lux 5 µm Amylose-2 250 x 4.6mm in NP

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 0.1% TFA in Hexane
Eluent B: 0.1% TFA in IPA
Gradient

Step No.	Time (min)	Pct A	Pct B
1	20	80	20

Profile:

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.
- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470
- Holder Part No.: KJ0-4282

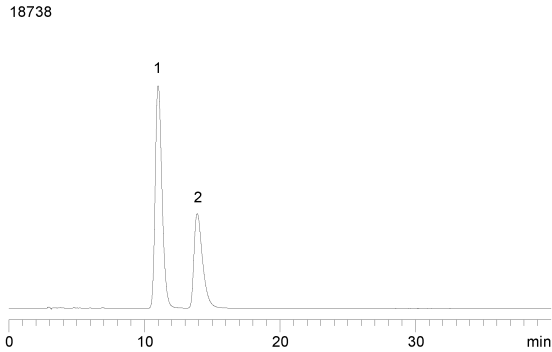

ANALYTES:

- 1 FMOC-Pro-OH
Retention Time: 11 min
- 2 FMOC-Pro-OH
Retention Time: 13.8 min

