

FMOC-DL-Leucine-(boc)-OH on Lux 5 μ m Amylose-2 250 x 4.6mm in RP

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: 0.1% TFA in CH₃CN

Eluent B: 0.1% TFA in Water

Gradient

Step No.	Time (min)	Pct A	Pct B
----------	------------	-------	-------

Profile:

1	10	60	40
---	----	----	----

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

18743

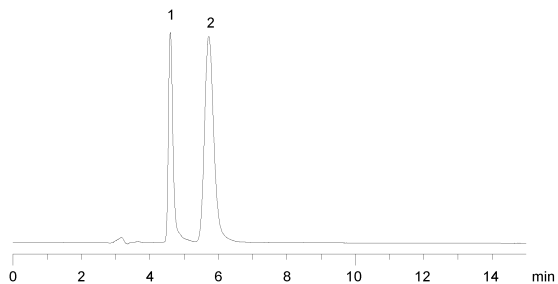

ANALYTES:

1 FMOC-Leu-OH

Retention Time: 4.5 min

2 FMOC-Leu-OH

Retention Time: 5.7 min

