

HPLC Application

ID No.: 18971

Propafenone on Lux 5 µm Amylose-2 by Chiral SFC

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Ethanol

Eluent B: Carbon Dioxide

Gradient

Step No.	Time (min)	Pct A	Pct B
----------	------------	-------	-------

Profile:

1	0	25	75
---	---	----	----

Flow Rate: 2.5 mL/min

Col. Temp.: 35 °C


Detection: UV-Vis Abs.-Diode Array (PDA) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

18971

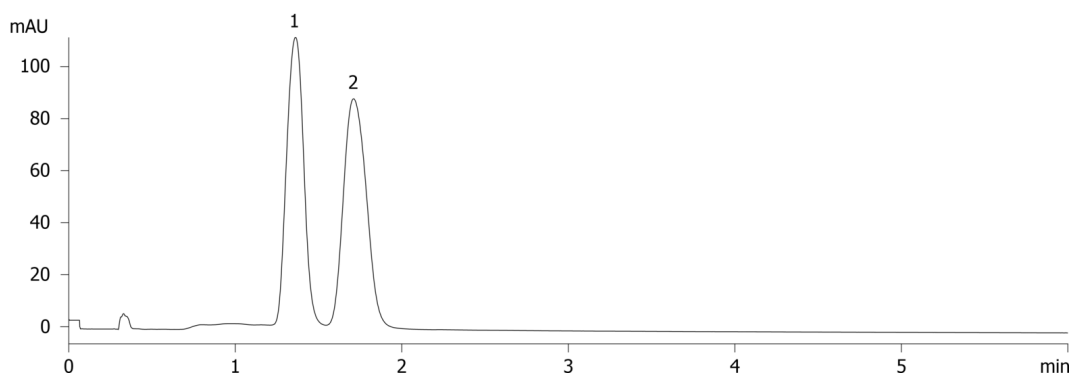

ANALYTES:

1 Propafenone

2 Propafenone

