

GC Application

ID No.: 19260

PolyWAX 850 on ZB-1XT 5x0.53x0.15

Column: Zebron[™] ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.15 μ m, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 5 meters x 0.53 mm x 0.15 μ m

Order No: 7AK-G026-05

Oven Profile: 35 °C to 450 °C @ 15 °C/min for 10 min

Carrier Gas: Constant Flow Helium, 15 mL/min

Injection: On-Column :1 1 μ L @ 38°C

Detection: Flame Ionization (FID) (450°C)


Analyst Note: Sample was 0.1% in CS₂. Chromatogram is baseline subtracted.



Products used in this application:



19260

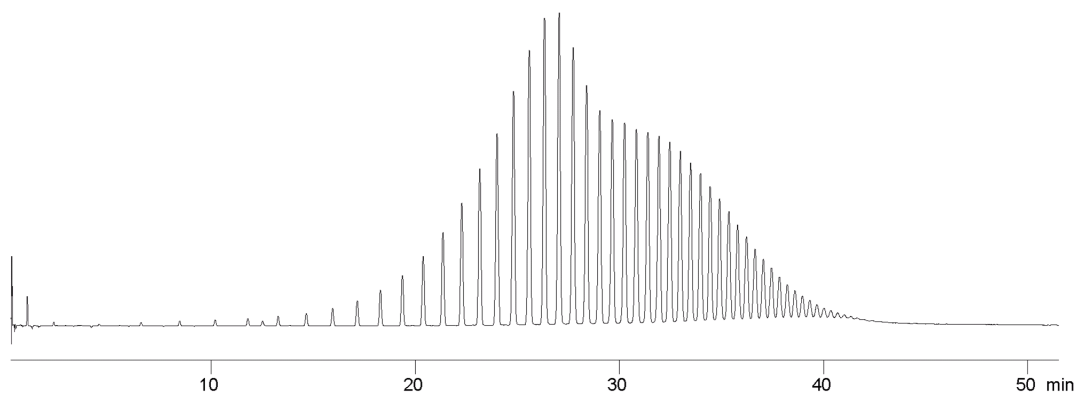

ANALYTES:

- 1 PolyWAX 850

