

ZB-1XT SimDist Test Mix

Column: Zebtron[™] ZB-1XT SimDist Metal, GC Cap. Column 5 m x 0.53 mm x 0.15 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 5 meters x 0.53 mm x 0.15 µm

Order No: 7AK-G026-05

Oven Profile: 105 °C Isothermal

Carrier Gas: Constant Pressure Helium, 5 mL/min

Injection: Split 20:1 0.05 µL @ 250°C

Detection: Flame Ionization (FID) (325°C)


Analyst Note: Sample at 1000µg/mL each in Hexane



Products used in this application:



19305

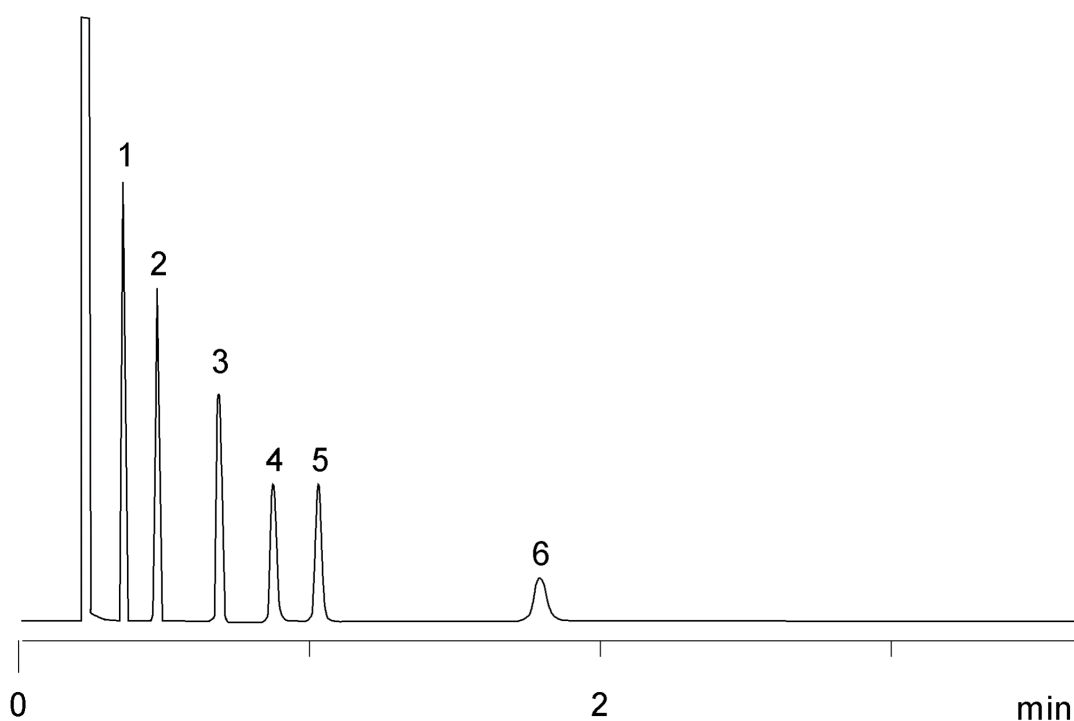

ANALYTES:

- 1 Undecane
- 2 Dodecane
- 3 Tridecane
- 4 1-Undecanol
- 5 Dicyclohexylamine
- 6 Pentadecane

