

ABN16 (acids, bases, neutrals) Mix using Kinetex 2.6µm C8, 50x2.1mm

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4497-AN

Elution Type: Gradient

Eluent A: 0.1% formic acid/water

Eluent B: 0.1% formic acid/ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.2	95	5
	3	4.2	5	95
	4	4.21	95	5
	5	5.5	95	5

Flow Rate: 0.8 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C8 for 2.1mm ID Columns, 3/Pk Part No.: AJ0-8784

- Holder Part No.: AJ0-9000



Products used in this application:

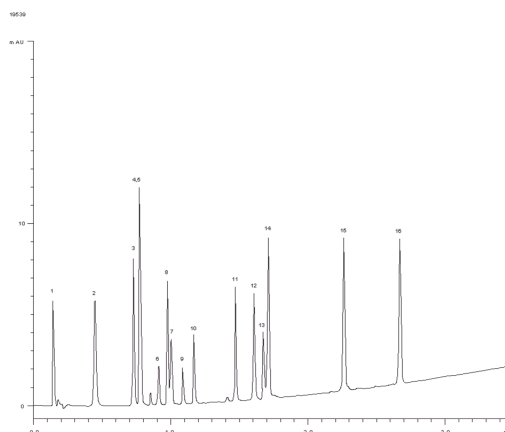

ABN16 (acids, bases, neutrals) Mix using Kinetex 2.6 μ m C8, 50x2.1mm

ANALYTES:

- 1 Pyridine
- 2 Acetaminophen
- 3 Sulfathiazole
- 4 Pindolol
- 5 Quinidine
- 6 Benzyl alcohol (Benzenemethanol)
- 7 Acebutolol
- 8 Phenol
- 9 Chlorpheniramine
- 10 Triprolidine (Actidil)
- 11 Prednisolone
- 12 3-Methyl-4-nitrobenzoic acid
- 13 Nortriptyline
- 14 2-Hydroxy-5-methylbenzaldehyde
- 15 Diflunisal
- 16 Hexanophenone

