

USP Norethindrone, Progesterone and Mestranol on Kinetex C8 2.6u 75x4.60mm @1.000mL/min

Column: Kinetex® 2.6 µm C8 100 Å, LC Column 75 x 4.6 mm, Ea

Dimensions: 75 x 4.6 mm ID

Order No: 00C-4497-E0

Elution Type: Isocratic

Eluent A: 50:50 Acetonitrile:Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	15	100

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C8 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8770
- Holder Part No.: AJ0-9000



Products used in this application:

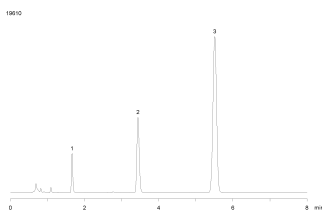

ANALYTES:

- 1** Norethindrone
Retention Time: 1.632 min
- 2** Progesterone
Retention Time: 3.439 min
- 3** Mestranol
Retention Time: 5.511 min

