

# HPLC Application

ID No.: 19818

## BSA Tryptic Digest on Aeris WIDEPORÉ 3.6u XB-C18

**Column:** Aeris<sup>™</sup> 3.6  $\mu$ m WIDEPORÉ XB-C18 200 Å, LC Column 150 x 2.1 mm, Ea

**Dimensions:** 150 x 2.1 mm ID

**Order No:** 00F-4482-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA in Water

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35
	2	3	65	35
	3	33	35	65

**Flow Rate:** 0.3 mL/min

**Col. Temp.:** 40 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)


**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORÉ C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000



Products used in this application:

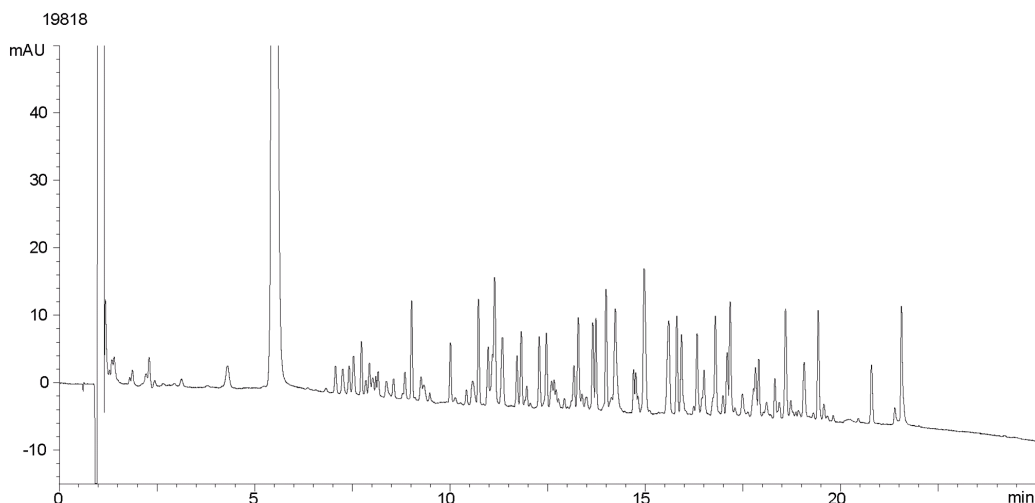

### ANALYTES:

- 1 Bovine Serum Albumin (BSA) (tryptic digest)

