

## Column Guard Competitive Test: Acquity Guard on Kinetex 1.7um XB-C18 50x2.1mm ID

**Column:** Kinetex® 1.7 µm XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea

**Dimensions:** 50 x 2.1 mm ID

**Order No:** 00B-4498-AN

**Elution Type:** Isocratic

**Eluent A:** ACN

**Eluent B:** Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

**Flow Rate:** 0.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)

**Analyst Note:** SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C18 for 2.1mm ID Columns, 3/Pk Part No.: AJ0-8782


- Holder Part No.: AJ0-9000



Products used in this application:



19975

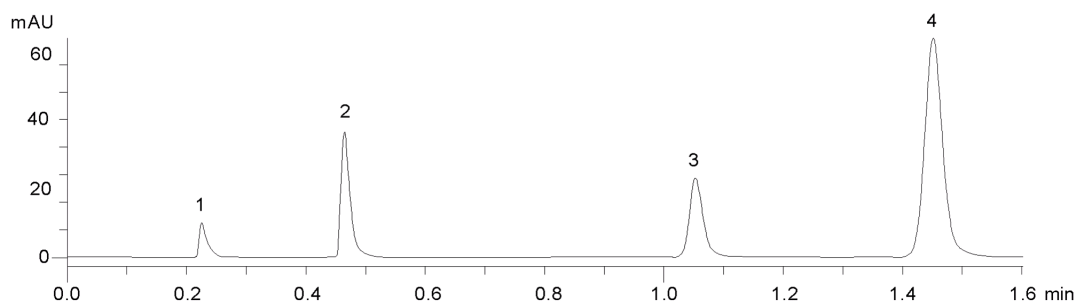

### ANALYTES:

- 1 Acetophenone
- 2 Benzene
- 3 Toluene
- 4 Naphthalene

