

Steroids Mix on Kinetex 2.6µm Phenyl-Hexyl

Column: Kinetex® 2.6 µm Phenyl-Hexyl 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4495-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	6	40	60

Flow Rate: 1.85 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC Phenyl for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8774

- Holder Part No.: AJ0-9000

19991

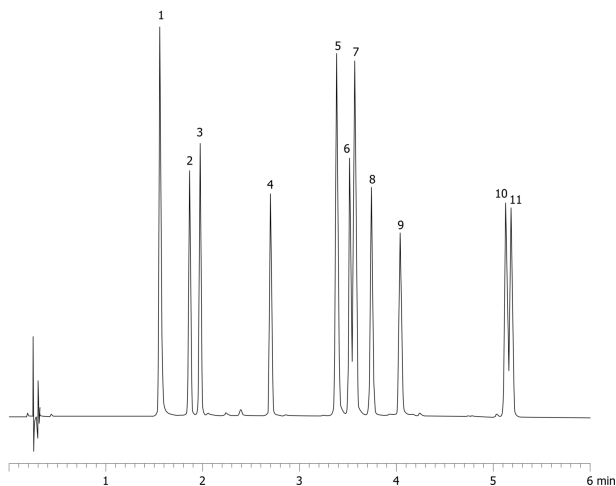

Products used in this application:



Steroids Mix on Kinetex 2.6 μ Phenyl-Hexyl

ANALYTES:

- 1 Estriol
- 2 Hydrocortisone
- 3 Prednisone
- 4 Corticosterone
- 5 beta-Estradiol
- 6 Cortisone acetate
- 7 alpha-Estradiol
- 8 21-Hydroxyprogesterone
- 9 Estrone
- 10 Progesterone
- 11 Deoxycorticosterone

