

## Ambucetamid on Lux 5µm Amylose-2 in NP

**Column:** N/A

**Dimensions:** N/A

**Order No:** N/A

**Elution Type:** N/A

**Eluent A:** Hexane / IPA (80:20) DEA (0.1%)

**Gradient**

Step No.	Time (min)	Pct A
----------	------------	-------

**Profile:**

1	0	100
---	---	-----

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

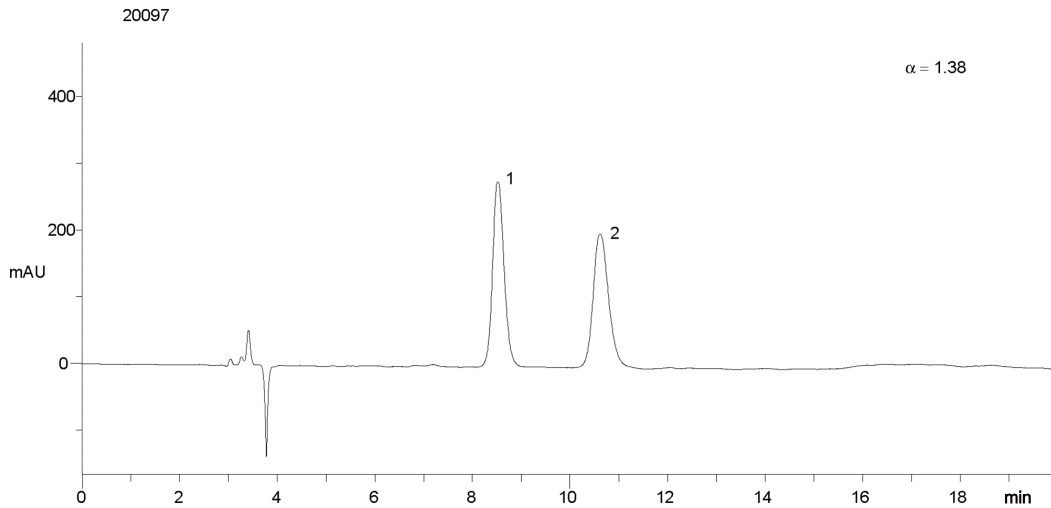

### ANALYTES:

**1** Ambucetamid

Retention Time: 8.5 min

**2** Ambucetamid

Retention Time: 10.6 min

