

# HPLC Application

ID No.: 20304

## Metomidate on Lux 5µm Cellulose-3 in RP

**Column:** Lux<sup>®</sup> 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4493-E0

**Elution Type:** Isocratic

**Eluent A:** MeOH with +0.1%DEA

**Eluent B:** 20mM NH<sub>4</sub>HCO<sub>3</sub>

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622

- Holder Part No.: KJ0-4282

20304



Products used in this application:

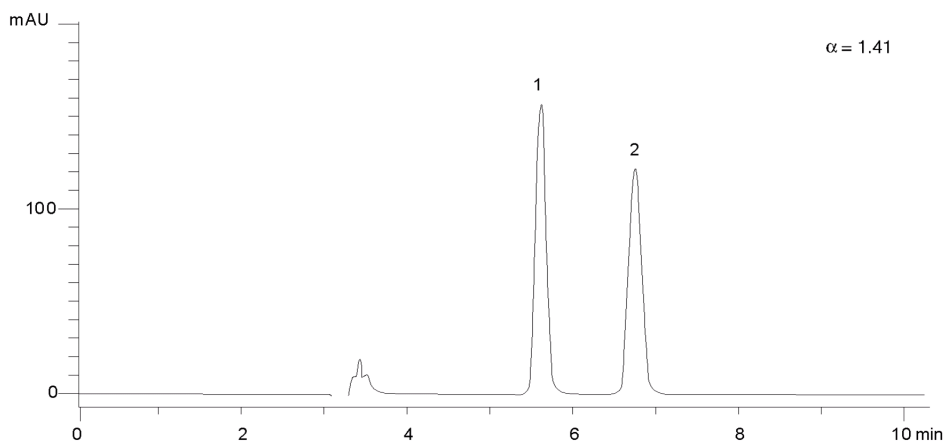

### ANALYTES:

- 1** Metomidate  
Retention Time: 5.62 min
- 2** Metomidate  
Retention Time: 6.75 min