

## Cobalt(III) acetylacetonate on Lux 5µm Amylose-2 in NP


**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/ EtOH (98:2) DEA (0.1%)**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

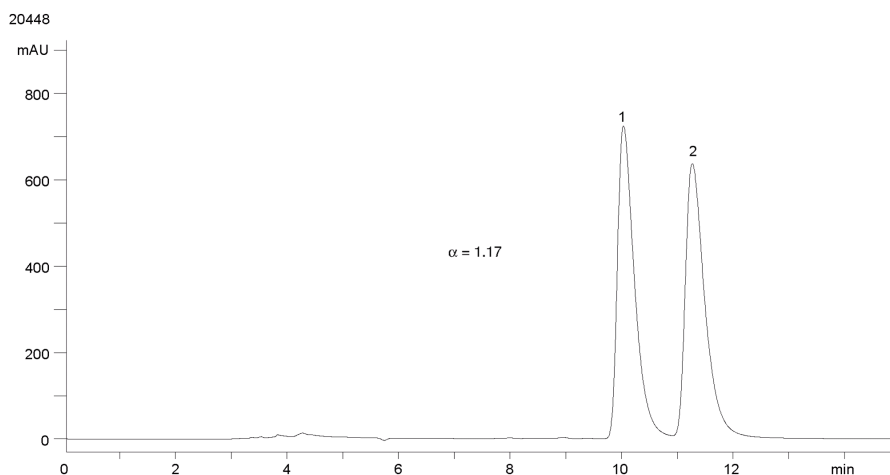

### ANALYTES:

**1** Cobalt (III) acetylacetonate

Retention Time: 10.04 min

**2** Cobalt (III) acetylacetonate

Retention Time: 11.27 min

