

HPLC Application

ID No.: 20451

Atropine on Lux 5 μ m Amylose-2 in NP

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane / EtOH (80:20) DEA (0.1%)

Gradient

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
|----------|------------|-------|

Profile:

| | | |
|----------|---|-----|
| 1 | 0 | 100 |
|----------|---|-----|

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

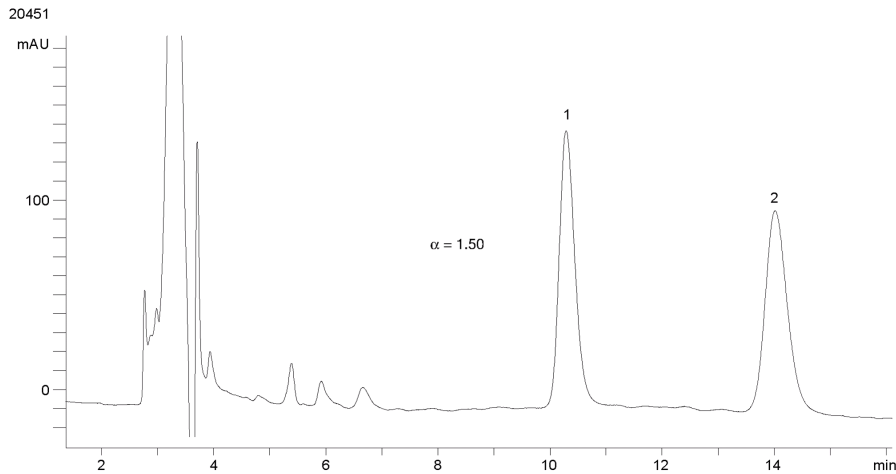

ANALYTES:

1 Atropine

Retention Time: 10.29 min

2 Atropine

Retention Time: 14.01 min

